

---

**BioAmp v1.5**

**Upside Down Labs**

**Jul 14, 2025**

# CONTENTS

|   |           |
|---|-----------|
| <b>1 BioAmp v1.5</b>  | <b>1</b>  |
| 1.1 Overview . . . . .  | 1         |
| 1.2 Features & Specifications . . . . .                       | 2         |
| 1.3 Hardware . . . . .  | 2         |
| 1.4 Contents of the kit . . . . .                             | 3         |
| 1.5 Software requirements . . . . .                           | 3         |
| 1.6 Assembling the kit . . . . .                              | 3         |
| 1.7 Using the kit . . . . .                                   | 6         |
| 1.7.1 Step 1: Connecting the cables . . . . .                 | 7         |
| 1.7.2 Step 2: Skin Preparation . . . . .                      | 7         |
| 1.7.3 Step 3: Electrode placements . . . . .                  | 7         |
| 1.7.4 Step 4: Connecting 9V battery . . . . .                 | 8         |
| 1.7.5 Step 5: Listen to your muscle signals . . . . .         | 9         |
| 1.7.6 Step 6: Visualize EMG signals on mobile phone . . . . . | 10        |
| 1.7.7 Step 7: Visualize the EMG signals on laptop . . . . .   | 12        |
| <b>2 Skin Preparation Guide</b>                               | <b>14</b> |
| 2.1 Why skin preparation is important? . . . . .              | 14        |
| 2.2 Kit Contents . . . . .                                    | 14        |
| 2.3 Steps to follow . . . . .                                 | 15        |
| 2.3.1 Step 1: Identify the targeted area . . . . .            | 15        |
| 2.3.2 Step 2: Apply NuPrep gel . . . . .                      | 15        |
| 2.3.3 Step 3: Clean the skin surface . . . . .                | 15        |
| 2.3.4 Step 4: Wipe off the gel . . . . .                      | 19        |
| 2.3.5 Step 5: Measuring the signals . . . . .                 | 19        |
| <b>3 Using BioAmp Bands</b>                                   | <b>22</b> |
| 3.1 Overview . . . . .  | 22        |
| 3.2 Why use BioAmp Bands? . . . . .                           | 22        |
| 3.3 Types of BioAmp Bands . . . . .                           | 22        |
| 3.3.1 1. Muscle BioAmp Band . . . . .                         | 23        |
| 3.3.2 2. Heart BioAmp Band . . . . .                          | 23        |
| 3.3.3 3. Brain BioAmp Band . . . . .                          | 24        |
| 3.4 Using Muscle BioAmp Band . . . . .                        | 25        |
| 3.4.1 Assembly . . . . .                                      | 25        |
| 3.4.2 Skin Preparation . . . . .                              | 27        |
| 3.4.3 Measure EMG . . . . .                                   | 27        |
| 3.5 Using Heart BioAmp Band . . . . .                         | 29        |
| 3.5.1 Skin Preparation . . . . .                              | 29        |
| 3.5.2 Assembly . . . . .                                      | 29        |
| 3.5.3 Measure ECG . . . . .                                   | 31        |
| 3.6 Using Brain BioAmp Band . . . . .                         | 32        |
| 3.6.1 Assembly . . . . .                                      | 32        |
| 3.6.2 Skin Preparation . . . . .                              | 32        |

|                                       |    |
|---------------------------------------|----|
| 3.6.3 Measure 1-channel EEG . . . . . | 32 |
|---------------------------------------|----|

---

CHAPTER  
ONE

---

## BIOAMP V1.5

2023 Edition

### 1.1 Overview

It is a small size portable biopotential amplifier with a no code setup to record and listen to your muscle signals (EMG) non invasively. The best part is that it doesn't require any microcontroller (like Arduino) to sample the signal. You just plug a 9V battery into the board, electrodes to the body, and an audio jack to a mobile/laptop, and you are ready to record signals from muscles (EMG) using softwares like audacity or Backyard Brain's spike recorder app.



## 1.2 Features & Specifications

|                       |                                       |
|-----------------------|---------------------------------------|
| Minimum Input Voltage | 7-9 V                                 |
| Input Impedance       | $10^7$ ohm                            |
| Fixed Gain            | $\sim 200x$                           |
| BioPotentials         | EMG (Electromyography)                |
| No. of channels       | 1                                     |
| Electrodes            | 3 (Positive, Negative, and Reference) |
| Dimensions            | 5.0 x 3.0 cm                          |
| Open Source           | Hardware + Software                   |

## 1.3 Hardware

Images below shows a quick overview of the hardware design.

### PCB Front



### PCB Back

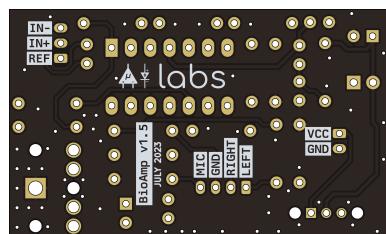


Fig. 1: Assembled BioAmp v1.5

## 1.4 Contents of the kit



BioAmp v1.5 bare PCB



BioAmp Cable v3



Muscle BioAmp Band



Components



9V Snap Cable



Gel electrodes



BioAmp AUX Cable

## 1.5 Software requirements

Before you start using the kit, download Backyard Brains' [Spike Recorder](#) or [Audacity](#) according to the operating system you are using (Windows, OSX, Linux).

## 1.6 Assembling the kit

You can get your own BioAmp v1.5 bag of parts from our [online stores](#) (shipping worldwide) and for assembling the board you can take a look at the step by step guide given below.

**Note:** Follow the highlighted yellow shapes to assemble your BioAmp v1.5!

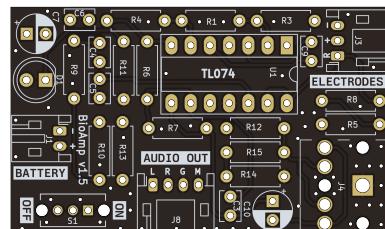


Fig. 2: Step 1 - Bare Board



Fig. 3: Step 2 - 100K Resistors



Fig. 4: Step 3 - 2.2K Resistor



Fig. 5: Step 4 - 1K Resistors



Fig. 6: Step 5 - 220K Resistors



Fig. 7: Step 6 - 10K Resistors



Fig. 8: Step 7 - 100nF Capacitors



Fig. 9: Step 8 - 1nF Capacitors



Fig. 10: Step 9 - JST PH Connectors



Fig. 11: Step 10 - IC Socket



Fig. 12: Step 11 - IC



Fig. 13: Step 12 - Power LED



Fig. 14: Step 13 - 47uF Capacitors

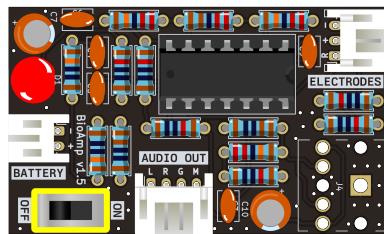


Fig. 15: Step 14 - Switch


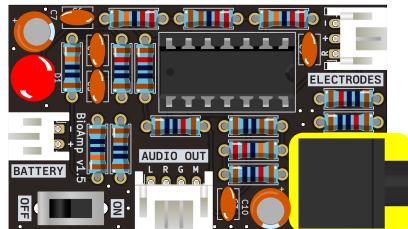


Fig. 16: Step 15 - Headphone jack


## 1.7 Using the kit

The image below shows the possibilities of using BioAmp v1.5. Seems complicated? Don't worry, we'll explain each and every step in detail. So follow along to create your own setup.

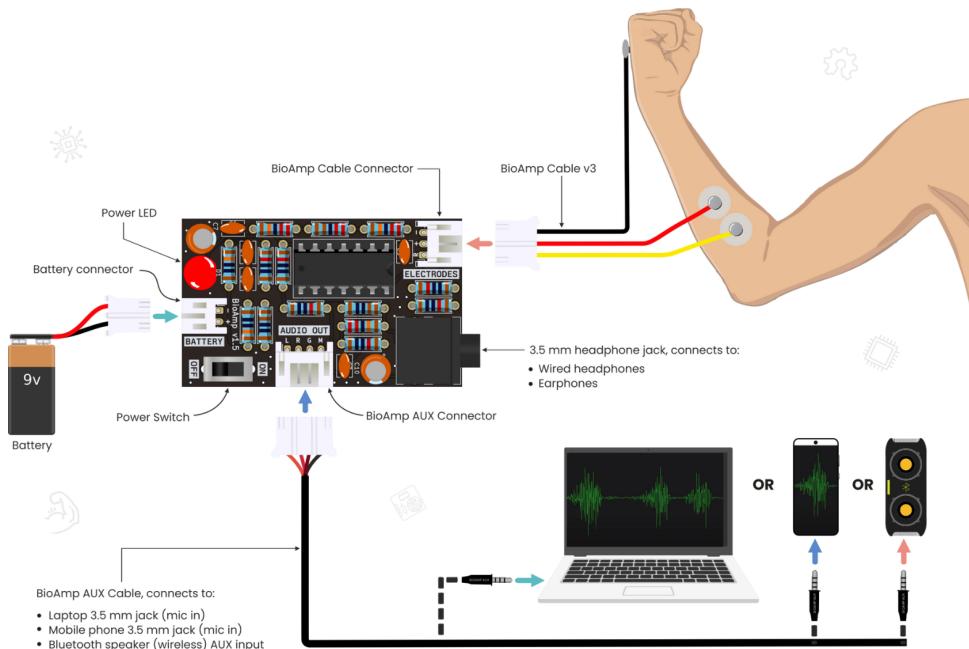

### 1.7.1 Step 1: Connecting the cables

Connect the BioAmp cable, 9V snap cable and BioAmp AUX cable to BioAmp v1.5 by inserting the cable ends in the respective JST PH connectors as shown below.

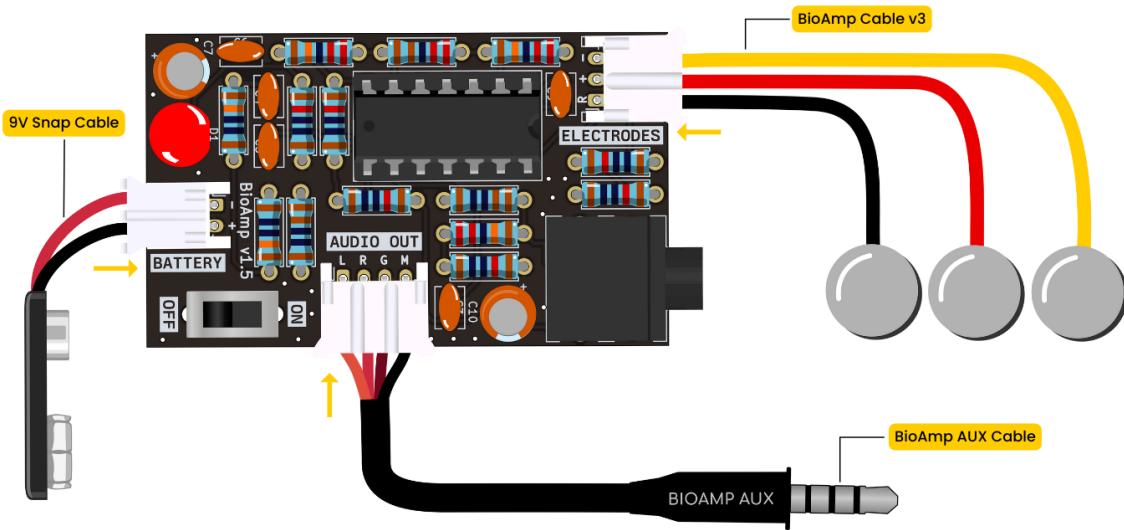

### 1.7.2 Step 2: Skin Preparation

Apply Nuprep Skin Preparation Gel on the skin surface where electrodes would be placed to remove dead skin cells and clean the skin from dirt. After rubbing the skin surface thoroughly, clean it with an alcohol wipe or a wet wipe.

For more information, please check out detailed step by step [Skin Preparation Guide](#).

### 1.7.3 Step 3: Electrode placements

We have 2 options to measure the EMG signals, either using the gel electrodes or using dry electrode based Muscle BioAmp Band. You can try both of them one by one.

#### Using gel Electrodes

1. Connect the BioAmp cable to gel electrodes.
2. Peel the plastic backing from the electrodes.
3. Place the IN+ and IN- cables on the arm near the ulnar nerve & REF (reference) at the back of your hand as shown in the connection diagram.



### Using Muscle BioAmp Band

1. Connect the BioAmp cable to Muscle BioAmp Band in a way such that IN+ and IN- are placed on the arm near the ulnar nerve & REF (reference) on the far side of the band.
2. Now put a small drop of electrode gel between the skin and metallic part of BioAmp cable to get the best results.

#### Note

These BioAmp Cable connections to gel electrodes/band are specifically for this version of BioAmp v1.5 only. The connections for other BioAmp hardware can be different.

#### Tip

- In this demonstration we are recording EMG signals from the ulnar nerve, but you can record EMG from other areas as well (biceps, triceps, legs, jaw etc) as per your project requirements. Just make sure to place the IN+, IN- electrodes on the targeted muscle and REF on a bony part.
- You can visit the complete documentation on how to [assemble and use the BioAmp Bands](#).

### 1.7.4 Step 4: Connecting 9V battery

Connect any 9V battery to BioAmp v1.5 using the 9V snap cable. Now activate the board by flipping ON the power switch, and you'll notice an LED light up, showing that the board is ready to use.



### 1.7.5 Step 5: Listen to your muscle signals

You can either listen to the muscle signals (EMG) on a speaker or wired earphones/headphones. Let's try both of them.

#### Using speakers

1. Connect the BioAmp AUX cable on a bluetooth speaker that have 3.5mm jack support.
2. Switch on the speaker and turn the volume to maximum.
3. Flex and listen to your muscles.



**Using wired earphones/headphones**

1. Plug your wired earphones or headphones on the 3.5mm jack of BioAmp v1.5.
2. Plug it in your ears.
3. Flex and listen to your muscles.

**1.7.6 Step 6: Visualize EMG signals on mobile phone**

Connect the BioAmp AUX cable to your mobile phone/tablet that has 3.5mm jack support. Now you can use various apps to visualise the signals.

**Using Phone Recorder app**

1. Open any audio recorder app on your mobile tablet.
2. Flex your muscle to be able to record the muscle signals.
3. If you want to extract that data then it will be saved by default as a .wav file but you can convert it in any other format according to your project requirements.



### Using Backyard Brains' Spike Recorder mobile app

1. Download the [Spike Recorder App](#) from playstore.
2. Open the app, click the setting icon on the top right corner and set the recording type to EMG.
3. Apply the 50Hz or 60Hz notch filter depending on the country you are living in. For example if you are in India then the AC current oscillates at a frequency of 50Hz but it oscillates at 60Hz frequency in USA. This AC current acts as a noise in the signals so we have to remove it by applying this notch filter.



4. Again click on the setting icon to close it and you are ready.
5. Flex your muscles to be able to visualize the muscle signals (EMG).
6. You can record the EMG data as a .wav file by pressing the record button on the top right corner of the app and then convert it in any other format as per your project requirements.



### 1.7.7 Step 7: Visualize the EMG signals on laptop

Connect the BioAmp AUX cable to your laptop that has 3.5mm jack support. Now you can install various softwares on your laptop to visualise the signals.

#### Using Backyard Brains' Spike Recorder

1. Install the spike recorder software that you downloaded earlier.
2. Open the software, click the setting icon on the top right corner and set the low band pass filter to 72Hz and high band pass filter to 720Hz.
3. Apply the 50Hz or 60Hz notch filter depending on the country you are living in. For example if you are in India then the AC current oscillates at a frequency of 50Hz but it oscillates at 60Hz frequency in USA. This AC current acts as a noise in the signals so we have to remove it by applying this notch filter.



4. Again click on the setting icon to close it and you are ready.
5. Flex your muscles to be able to visualize the muscle signals (EMG)
6. You can record the EMG data as a .wav file by pressing the record button on the top right corner of the app and then convert it in any other format as per your project requirements.



### Using Audacity

1. Install the audacity software that you downloaded earlier.
2. Open the software and set the recording device to microphone.
3. Flex your muscles to be able to visualize the muscle signals (EMG)
4. By default the EMG data would be recorded as a .wav file but you can convert it in any other format as per your project requirements.



---

**CHAPTER  
TWO**

---

## **SKIN PREPARATION GUIDE**

### **2.1 Why skin preparation is important?**

Proper skin preparation is crucial before recording any biopotential signal be it Electrocardiography (ECG), Electromyography (EMG), Electroencephalography (EEG), or Electrooculography (EOG).

- **Clean skin surface:** Removes dead skin cells, oils, & other substances that increases skin impedance.
- **Improve impedance:** Improves the conduction of electrical signals from the body to the recording equipment and minimizes impedance.
- **Electrode-skin contact:** Ensures optimal contact between the electrodes and the skin surface.
- **Signal quality:** Enhances the overall quality of recorded signals, providing clear & reliable data for analysis & improves the ability to capture subtle variations in biopotential signals.
- **Consistency in recordings:** Reduces variability in signal quality, making it easier to make any Human-Computer Interface (HCI), Brain-Computer Interface (BCI) project or a real-world application.
- **Long term adhesion:** Facilitates long-term adhesion & stable placement of electrodes to the skin during extended signal monitoring.

### **2.2 Kit Contents**

|                         |  |
|-------------------------|--|
| Nuprep gel              | Mildly abrasive, highly conductive gel that should be applied before placing the electrodes on the skin to improve signal quality & enhances the performance of monitoring electrodes. |
| Electrode Gel           | Highly conductive gel that acts as a coupling agent between dry electrodes and the skin to aid the transmission of biopotential signals like ECG, EMG, EOG, or EEG.                    |
| Ten20 paste             | Contains the right balance of adhesiveness and conductivity, enabling the dry electrodes to remain in place while allowing the transmittance of biopotential signals.                  |
| Alcohol Swabs/Wet wipes | Soft & non-woven pads that helps in cleaning the skin surface and does not leave any residue.  |
| Cotton Swabs            | Useful while applying nuprep gel or ten20 paste.   |

# Contents of the kit



NuPrep Gel



Cotton Swabs



Electrode Gel



Alcohol Swabs



Ten20 Paste

## 2.3 Steps to follow

You can follow the steps given below to do the skin preparation properly:

### 2.3.1 Step 1: Identify the targeted area

Identify the target area where the gel electrodes or BioAmp Bands will be placed for recording the biopotential signals.

### 2.3.2 Step 2: Apply NuPrep gel

Take a small amount of NuPrep gel using a cotton swab and apply it on your targeted area.

### 2.3.3 Step 3: Clean the skin surface

Use gentle, circular motions to rub the gel on the skin surface. This removes all the dead skin cells & improves conductivity.

#### Warning

Do not rub the gel for too long as it has abrasive properties and may cause skin redness and irritation.



Fig. 1: Target area to record EOG



Fig. 2: Target area to record EMG



Fig. 3: Target area to record ECG

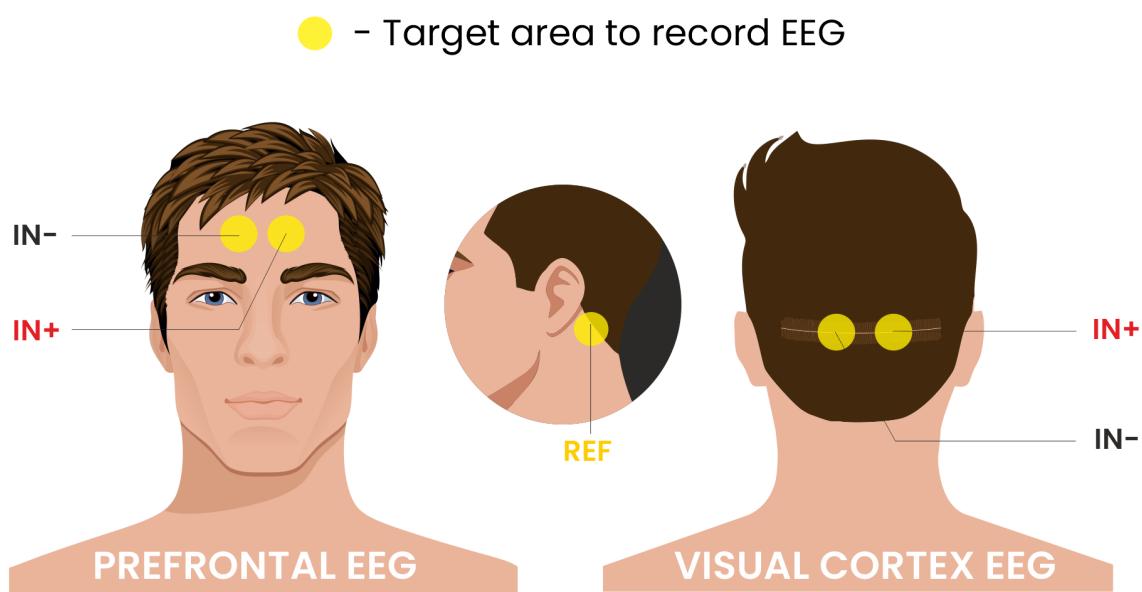


Fig. 4: Target area to record EEG



Fig. 5: Rub the gel gently using the cotton swab

### 2.3.4 Step 4: Wipe off the gel

Wipe away excess gel with alcohol swabs or wet wipes.



Fig. 6: Wipe away access gel

#### Warning

- Using alcohol swabs can dry out the skin, so don't use them if your skin is already dry.
- Close your eyes while using the alcohol swabs for EOG recording else it may cause eye redness & irritation.

### 2.3.5 Step 5: Measuring the signals

Now you can either use gel electrodes or BioAmp bands for the signal recording.

#### Using gel electrodes

Connect the BioAmp cable to gel electrodes, peel the plastic backing from electrodes and place the IN+, IN-, REF cables according to your specific biopotential recording.

#### Note

While placing the gel electrodes on the skin, make sure to place the non-sticky tab of the electrode in the direction opposite to your hair growth. This allows you to remove the electrodes easily without pulling off much body hair.



Fig. 7: Placing gel electrodes on skin surface

### Using BioAmp bands

Connect the BioAmp cable to your BioAmp band. Now apply a small amount of **electrode gel** or **Ten20 conductive paste** on the dry electrodes between the skin and metallic part of BioAmp cable. This improves the signal conductivity, enhancing overall signal quality.

#### Note

The above graphics demonstrates the use of electrode gel/Ten20 paste with Muscle BioAmp Band. Similarly you can use Brain BioAmp Band and Heart BioAmp Band. Refer to [\*Using BioAmp Bands\*](#) guide to assemble and use all the BioAmp Bands correctly.

Now you are all set! Make all the connections correctly and start recording your biopotential signals.

#### Warning

NuPrep gel, Ten20 paste and the alcohol swabs shouldn't be used if you have a history of skin allergies to lotions and cosmetics.



Fig. 8: Method 1: Using Electrode gel



Fig. 9: Method 2: Using Ten20 paste

## USING BIOAMP BANDS

### 3.1 Overview

BioAmp Bands are dry electrode-based stretchable bands that allows you to record biopotential signals from your body be it from brain (EEG), muscles (EMG) or heart (ECG). These bands can only be used with our BioAmp Hardware by making the connections using BioAmp Cable.

### 3.2 Why use BioAmp Bands?

Usually, people use gel electrodes to record biopotential signals from the skin surface. But, it has its own disadvantages. So we came up with these BioAmp Bands using which users can enjoy a more comfortable, cost-effective, and hassle-free experience while recording biopotential signals.

- **Comfort:** BioAmp Bands are generally more comfortable to wear than gel electrodes, especially for long-term recordings. They conform to the body's shape and avoid the sticky, sometimes irritating sensation of gel electrodes.
- **Reusability:** Unlike gel electrodes, which are often single-use and need to be replaced frequently, BioAmp Bands can be reused multiple times. This makes them more cost-effective and environmentally friendly.
- **Ease of Use:** These bands are easy to wear and adjust, reducing the hassle of setup and ensuring consistent placement.
- **Hygiene:** They can be easily cleaned and sanitized between uses, reducing the risk of skin irritation and infections. Gel electrodes, on the other hand, can leave residue on the skin surface.
- **Performance:** The bands can provide stable and reliable signal recordings depending on your environment conditions. For hot/humid conditions, the bands usually perform better while recording the signals. But if the weather is cold causing dry skin, then it is recommended to prepare the skin properly and apply electrode gel between the metallic part of cable and skin surface. If you feel that the skin impedance is increasing, then reapply electrode gel frequently. The other option is to use gel electrodes after preparing the skin properly.

### 3.3 Types of BioAmp Bands

There are 3 types of BioAmp Bands and all these bands offer targeted and efficient solutions for recording biopotential signals from the muscles, heart, and brain, making them versatile tools for a wide range of HCI/BCI applications.

### 3.3.1 1. Muscle BioAmp Band

Muscle BioAmp Band (EMG Band) is a stretchable band that can be connected to any of our Muscle BioAmp Hardware or any EXG sensor using a BioAmp Cable. It allows you to record your muscle signals hassle-free.



|                     |   |
|---------------------|---|
| Length              | 13 inches                                 |
| Stretchability      | 2X (Upto 26 inches)                       |
| Usability           | Reusable as it comes with washable fabric |
| Interface           | Snap electrodes                           |
| Compatible Hardware | Muscle BioAmp Hardware or any EXG sensor  |
| BioPotentials       | EMG                                       |
| No. of channels     | 1   |
| Wearable            | Yes                                       |

### 3.3.2 2. Heart BioAmp Band

Heart BioAmp Band (ECG Band) is a stretchable band that can be connected to any of our Heart BioAmp Hardware or any EXG sensor using BioAmp Cable. It allows you to record your ECG signals hassle-free.



|                     |   |
|---------------------|---|
| Length              | 37 inches                                 |
| Stretchability      | 2X (Upto 74 inches)                       |
| Usability           | Reusable as it comes with washable fabric |
| Interface           | Snap electrodes                           |
| Compatible Hardware | Heart BioAmp Hardware or any EXG sensor   |
| BioPotentials       | ECG                                       |
| No. of channels     | 1   |
| Wearable            | Yes                                       |

### 3.3.3 3. Brain BioAmp Band

Brain BioAmp Band (EEG Band) is a stretchable band that can be connected to any of our Brain BioAmp Hardware or any EXG sensor using BioAmp Cable to record signals from the brain hassle-free.

|                     |   |
|---------------------|---|
| Length              | 15.5 inches                               |
| Stretchability      | 2X (Upto 31 inches)                       |
| Usability           | Reusable as it comes with washable fabric |
| Interface           | Snap electrodes                           |
| Compatible Hardware | Brain BioAmp Hardware or any EXG sensor   |
| BioPotentials       | EEG                                       |
| No. of channels     | 2 or 6                                    |
| Wearable            | Yes                                       |

You can get either a 2-channel or a 6-channel Brain BioAmp Band according to your project or research requirements:

#### 2-Channel Brain BioAmp Band

It can be used to record EEG signals up to 2 channels either from the visual cortex (back of your head) or the prefrontal cortex part of brain.



## 6-Channel Brain BioAmp Band

It can be used to record EEG signals up to 2 channels either from the visual cortex (back of your head) or the prefrontal cortex part of brain.



## 3.4 Using Muscle BioAmp Band

### 3.4.1 Assembly

1. Take your Muscle BioAmp Band, hold the side of the band that has buckle on it and align the top part of the buckle with the flat surface of the snap.



1. Take the other end of the band and insert it in the buckle.
3. Your band is now ready to use. You can also adjust the size of the band according to your targeted muscle.

Take other end of the band and  
insert it in the buckle



Insert this end of the band into the buckle  
as shown and pull it to finish the assembly





### 3.4.2 Skin Preparation

Apply Nuprep Skin Preparation Gel on the skin surface where dry electrodes would be placed to remove dead skin cells and clean the skin from dirt. After rubbing the skin surface thoroughly, clean it with an alcohol wipe or a wet wipe.

For more information, please check out detailed step by step [Skin Preparation Guide](#).

### 3.4.3 Measure EMG

1. Flip the band and snap the dry electrodes of the BioAmp Cable on it as shown below.



1. Flip the band again and wear it on your arm in such a way that IN+ and IN- are placed on the arm near the ulnar nerve and REF (reference) on the far side of the band.

**Note**

Make sure the dry electrodes (shiny parts of the BioAmp Cable) are in direct contact with the skin.

3. Now put a small amount of electrode gel or Ten20 paste between the skin and dry electrodes to get the best signal acquisition.



**Note**

- After using the band, don't leave the gel residue on the dry electrodes longer than an hour as it may corrode them over a period of time.
- Wash the band with liquid soap and rinse it properly after every use. Use it again only when it is completely dry.

## 3.5 Using Heart BioAmp Band

### 3.5.1 Skin Preparation

Apply Nuprep Skin Preparation Gel on your chest where dry electrodes would be placed to remove dead skin cells and clean the skin from dirt. After rubbing the skin surface thoroughly, clean it with an alcohol wipe or a wet wipe.

For more information, please check out detailed step by step [Skin Preparation Guide](#).

### 3.5.2 Assembly

1. Take your Heart BioAmp Band and wrap the band around your chest in such a way that the pointy part of the snap touches your chest and the flat part is on the outer side.



2. Now insert the loose end of the band into the buckle and tighten it by pulling the strap.
3. Your band is now ready to use. You can also adjust the size of the band according to your chest size.



### 3.5.3 Measure ECG

1. Snap the IN- cable on the left most side of the band, IN+ cable in the middle, and REF cable on the right side as shown below.



#### Note

Make sure the dry electrodes (shiny parts of the BioAmp Cable) are in direct contact with the skin.

2. Now put a small amount of electrode gel or Ten20 paste between the skin and dry electrodes to get the best signal acquisition.



**Note**

- After using the band, don't leave the gel residue on the dry electrodes longer than an hour as it may corrode them over a period of time.
- Wash the band with liquid soap and rinse it properly after every use. Use it again only when it is completely dry.

## 3.6 Using Brain BioAmp Band

### 3.6.1 Assembly

You get the band in two parts - the longer part consists of buckles at both ends and the shorter one has loose ends on both sides.

1. Hold one end of the longer band and align the top part of the buckle with the flat surface of the snap.
2. Now take the shorter band and insert it into the buckle of longer band.
3. Repeat step 1 and 2 for the other buckle on the longer band.
4. Your band is now ready to use. You can also adjust the size of the band according to your head size.

### 3.6.2 Skin Preparation

Apply Nuprep Skin Preparation Gel on your targeted area (visual cortex or prefrontal cortex) where dry electrodes would be placed to remove dead skin cells and clean the skin from dirt. After rubbing the skin surface thoroughly, clean it with an alcohol wipe or a wet wipe.

For more information, please check out detailed step by step [Skin Preparation Guide](#).

### 3.6.3 Measure 1-channel EEG

1. Flip the band, take your BioAmp Cable, and snap the REF cable on a gel electrode. Now snap the IN- and IN+ cable on:
  - Fp1 and Fp2 positions for recording EEG from prefrontal cortex
  - O1 and O2 positions for recording EEG from visual cortex

**Note**

The electrode positions mentioned above are according to [International 10-20 system for recording EEG](#).

2. Flip the band again and wear it in a way so that the dry electrodes (shiny parts of the cable) are in contact with:
  - skin surface on the forehead (if recording from prefrontal cortex)
  - scalp surface on the back side of your head (if recording from visual cortex)
3. Peel off the plastic backing of the gel electrode and place it on the bony part behind your earlobe.

**Note**

While placing the gel electrodes on the skin, make sure to place the non-sticky tab of the electrode in the direction opposite to your hair growth. This allows you to remove the electrodes easily without pulling off much body hair.

4. Now put a small amount of electrode gel or Ten20 paste between the skin/scalp and dry electrodes to get the best signal acquisition.